



# Product Briefing

## **Chromacol 05-MTVC-96 WebSeal™ Plates**

500µL Vials with Sealing Mat


- A new addition to the range of WebSeal™ products.
- New 500µL vial pre-inserted into a mid-depth micro plate.
- Fast and secure sealing with the WebSeal™ WSM-I silicone/PTFE mat.
- Compatible with High Throughput HPLC auto samplers.



Chromacol introduces a new addition to the WebSeal™ range of 96-Well products. The 05-MTPVC-96 consists of 5 mid-depth polypropylene plates pre-inserted with 96 500µL glass sample vials. Supplied complete with 5 WSM-I sealing mats for fast and secure closure of all vials in less than 60 seconds.

### **05-MTPVC-96 WebSeal™ plates provide:-**

- Broad compatibility with HTS and Well-plate auto samplers.
- PTFE coated high purity silicone formulation gives chemically inert chamber.
- Precision tapered base gives high sample recovery
- Lower profile gives access to majority of modular well-plate samplers.
- Mat is compatible with existing deep-well MTPVC-96.
- Individual vials may be removed with a special cutting tool.

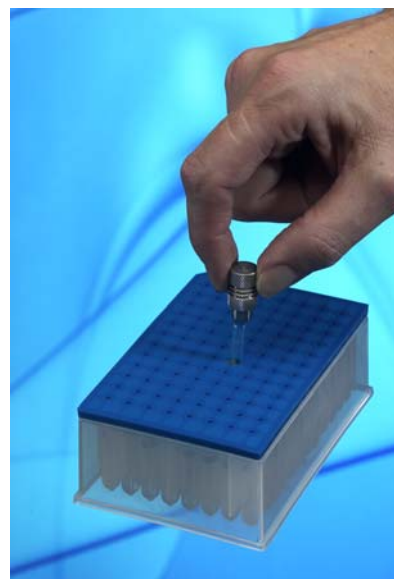

### **Chromacol 500uL WebSeal™ Products**

<b>Chromacol Part</b>	<b>Pack Size</b>	<b>Description</b>
05-MTPVC-96	5	500µL tapered vials pre-inserted into a 96 well mid-well plate with separate sealing mat.
05-MTVC-96	1000	500µL tapered vials 31mm x 6mm for 96-well plates
MTPC-I	1	WebSeal™ cutter for vial removal
WSM-I	5	96 well, silicone/PTFE coated elastomer WebSeal™ mat