



Product Briefing



Chromacol Autocrimp™

Pneumatic Crimpers.

Crimp-top vials are used extensively in automatic sampling systems - typically autosamplers for HPLC and GC.

Autocrimp™ makes the task of sealing large series of vials very easy due to the pneumatic operation of the crimper. Sealing a vial with a hand crimper is no problem with small numbers of vials, but is very hard on the wrist and lower arm when it comes to large numbers.

Autocrimp™ is a bench top model equipped with a foot pedal for activation and therefore leaves both hands free. The base plate is made of heavy steel for stability as is the vertical shaft. The base plate is also predrilled to allow permanent mounting.

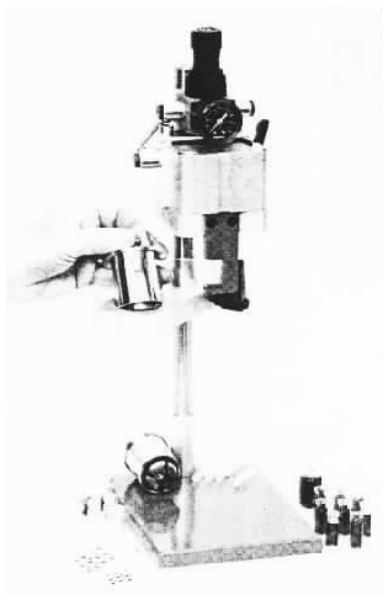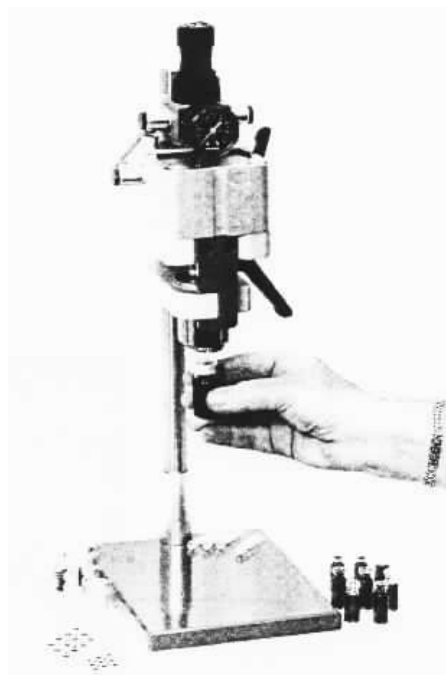

Autocrimp™ makes perfect seals with all types of crimp caps because the crimp jaws are exchangeable and adjustable.

Crimp jaws are available for 8, 11, 13 and 20 mm crimp caps and the stroke lengths can be fine tuned to obtain a perfect seal independent of the thickness of the sealing membrane.

For easy removal of 11 and 20mm crimp caps a special de-capping tool is available. It mounts in the same holder as the crimp jaw.

Specifications

Benchtop model with steel shaft and base, including foot pedal.

Pressure Required 5 - 6 Bar

Exchangeable 8, 11, 13 and 20mm crimping and de-capping jaws.

Height adjustable to minimum 55cm

Base Dim

Ordering

Chromacol Part Number	Description
CMSP-0	Pneumatic Benchtop Crimping Station
CMJ-8	8mm Crimp Jaws
CMJ-11	11mm Crimp Jaws
CMJ-13	13mm Crimp Jaws
CMJ-20	20mm Crimp Jaws